

Purite Integra L

Pure water systems for centralised distribution of laboratory water



Our Purite Integra L units have been designed to provide a broad spectrum of purified water for laboratory applications. Each unit is designed to feed a laboratory suite with several outlets using a ring main to provide a continuous supply. The water is stored in an integral 250 litre stainless steel tank. Two units are available: The Integra LH for hard feed water and Integra LS for soft feed water.

Additional features

- Single compact unit which purifies, stores and distributes lab water
- Produces up to 600 l/hr
- Can supply ring mains of potentially up to 200 metres (application and site dependent)
- Provides ASTM Type I, II or III water and BS EN ISO grade 1,2 or 3 water
- Utilises reverse osmosis technology, in conjunction with activated carbon and particulate filtration as part of the first stage of purification
- Integrated data logging for performance traceability equipped as standard
- Delivered factory tested for ease of installation
- LCD screen display for ease of operation.

Purite Integra L

Pure water systems for centralised distribution of laboratory water


Purite Integra L Features

Our Integra range is ideal for laboratories requiring greater daily volumes of purified water with several points of use.

Our Integra HP and Integra L systems deliver laboratory grade deionised water from 60-600 litres per hour, while our Integra 200E system is a low energy, self-contained unit utilising the latest low energy reverse osmosis membranes and electro-deionisation technology.

To provide higher grades of purified water the basic Integra L units can be supplemented by the addition of optional polishing packs comprising of stand-alone exchangeable cylinders. For applications requiring high quality water with enhanced bacterial specification, UV disinfection and 0.2µm filtration is available in an optional BIO pack.

Purite Integra L process flow

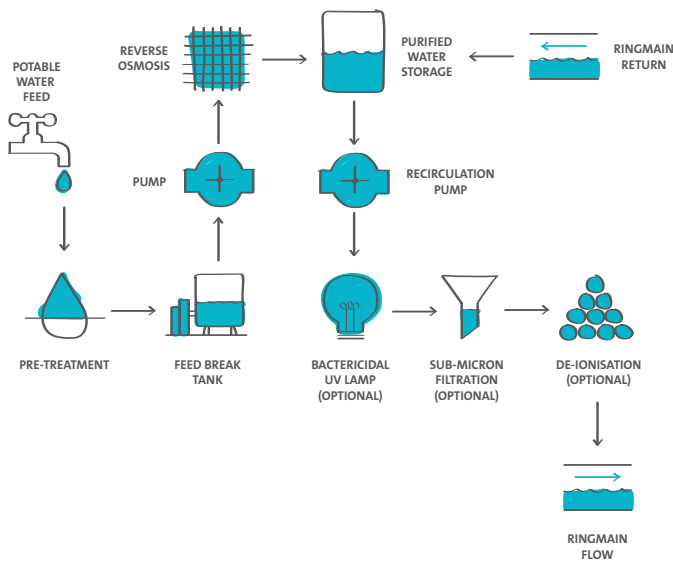

Visit online!

PURITE

Veolia Water Purification Systems Ltd trading as Purite

Bandet Way Thame Oxfordshire OX9 3SJ United Kingdom

T +44(0)1844 217141

E vtc.vwts.puritemail.all@veolia.com

W purite.com

Errors and omissions excluded. Purite reserves the right to change the specification in accordance with our programme of continual improvement.

Unit specifications

	Integra L
Width (mm)	1000
Depth (mm)	750
Height (mm)	1800
Max shipping weight (kg)	310
Max working weight (kg)	550
Installation requirements	
Power	Single Phase, 230V, +/- 10%, 50 Hz
Feed water	Potable*
Maximum TDS (ppm)	1000
Minimum inlet pressure - psi (bar)	15 (1)
Maximum inlet pressure - psi (bar)	90 (6.2)
Feedwater temperature	1 - 30°C
Feedwater flowrate (litres/hour)	-
Free chlorine	Must be dechlorinated

* Softened feed water required for 600 l/hr (Integra LS). Integra LH can operate on hard water up to 400ppm as CaCO₃

System specifications

	Integra L
Pure water storage	250 litres
Display panel - LCD colour touch screen	LCD screen
Pre-treatment	5µm pre-filter
Reverse osmosis	Low energy membranes
Deionisation	Optional Di packs*
Micro filtration	Optional 0.2µm**
UV lamp	Optional 254nm**
Purified water make-up flow rate @ 10°C	200 - 600 l/hr***
Purified water distribution	Up to 3000 l/hr
Carbon dioxide degassing	-
TOC reduction	Optional*

* 10, 15, 18MΩ.cm polishing deionisation packs available including activated carbon for TOC reduction

** Included in BioPack

*** Softened feed water required for 600 l/hr (Integra LS). Integra LH can operate on hard water up to 400ppm as CaCO₃

Treated water specifications

	Integra L
Conductivity	< 30µS/cm*
Bacteria	> 99% rejection**
Organics – TOC (ppb)	< 50***
Particles	0.2µm**
Endotoxins	0.25EU/ml**

* For enhanced inorganic and organic quality the unit can be fitted with either 10,15 or 18MΩ.cm polishing deionisers packs.

** Optional 'BioPack' will provide purified water with a total viable count of <1cfu/ml, endotoxin level <0.25EU/ml and particles <0.2µm

*** Typical level, depends on feedwater