


PURITE

# Purite Remote Dispenser

Ultrapure water at your fingertips

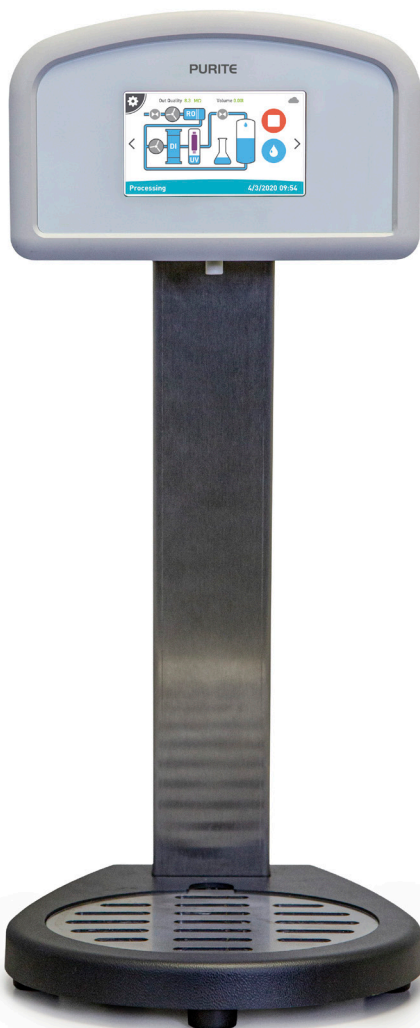

Discover seamless efficiency and unrivalled convenience with our Remote Dispenser, the latest addition to our Purite bench-top water purification systems.

Our state-of-the-art 'satellite' dispenser allows you to access high-quality water directly at your workstation, eliminating the need to reach out to the main RO system.

## Features:

- Space-saving - compact and practical design
- Ergonomic workflow - improved user comfort
- Flexibility - can be easily repositioned as your lab layout evolves
- Reach - Position up to 5 meters from the main water purification system
- Reliability - backed by Purite's reputation for quality and durability
- Seamless integration - effortlessly connect to your existing Purite system
- High quality water - matches water quality produced by the main Purite system
- Integrated UV LED - microbiological inactivation at the point of dispense
- LCD colour touch screen panel.

**WATER PURIFICATION SYSTEMS**


# Purite Remote Dispenser

Ultrapure water at your fingertips

In demanding laboratory environments, every movement matters and ergonomics play a crucial role in optimising productivity and comfort. Our extremely compact Remote Dispenser can be conveniently and discreetly positioned anywhere on your bench, saving valuable lab space.

The Purite Remote Dispenser is designed with multiple ergonomic and practical features ensuring total convenience and efficiency along with versatile installation, storage and dispensing options for any workstation.

## Purite Remote Dispenser process flow\*

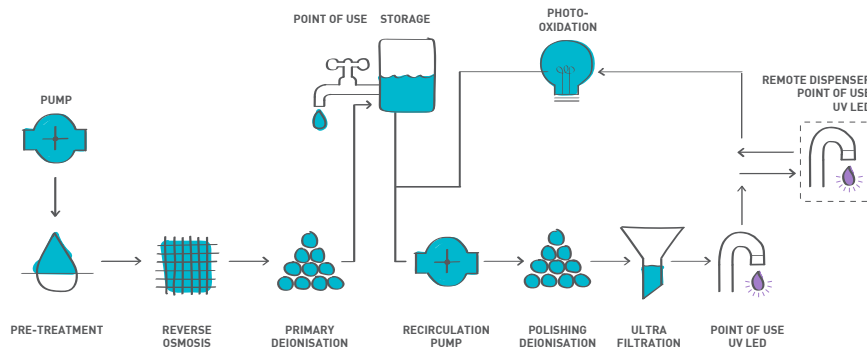

\* The process flow diagrams differ according to which water purification system the dispenser is connected to.

## Features & technical specifications

Unit features & technical specifications			
Model No.	L300940		
Dimensions	Height: 524 mm	Width: 228 mm	Depth: 305 mm
Removable drip tray	Yes		
Point of use mercury free UV LED water disinfection	Yes		
HMI interactions	Mirror's screen functionality on the main water purification system		
Water quality	Matches quality produced from main water purification system		
Reachable distance from feed water system	3 metre and 5 metre connections available		
Number of remote dispensers operated from one Purite lab water purification system	1		
Simultaneous product water dispensing	No		
Remote dispenser compatibility	Purite Fusion WiFi, Purite HP		
Connections provided by main water purification system	Feed water in - 8mm tube, Feed water return - 8mm tube, USB & POD connection		
Electrical requirements	Electrical supply: 24V   Max power consumption: 1 Amp (Provided by main system)		
Materials	Polypropylene base. Stainless steel 316 stand & removable grill. Polypropylene HMI pod		

## PURITE

Veolia Water Purification Systems Ltd trading as Purite

Bandet Way Thame Oxfordshire OX9 3SJ United Kingdom

T +44 (0)1844 217141

E vtc.vwts.puritemail.all@veolia.com

W purite.com

Visit online!

